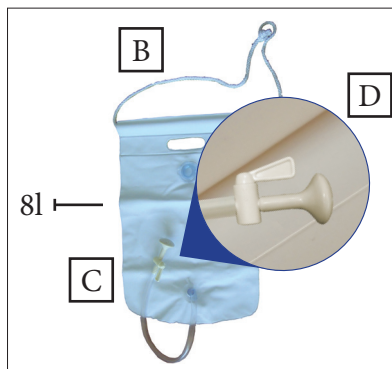
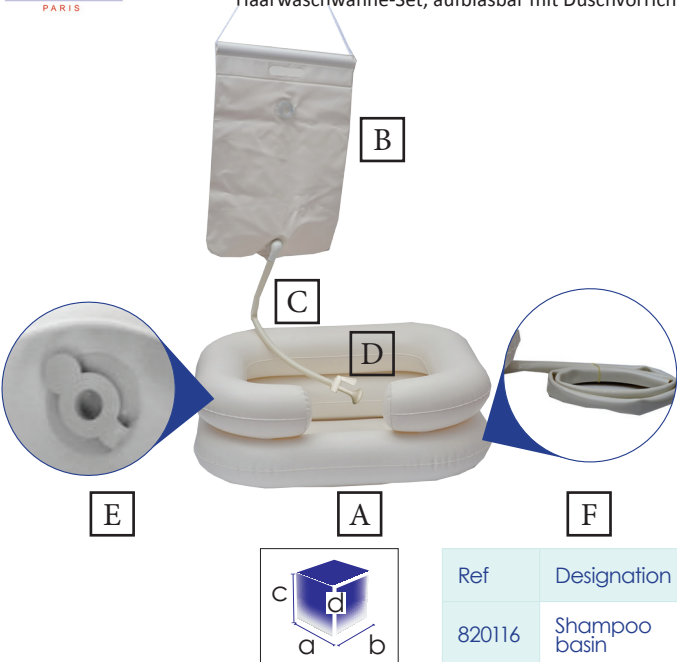
Ref	Designation	a	b	c	d	Weight
820116	Shampoo basin	61 cm	53 cm	20 cm	10l	0.83 kg

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EN- Please read these safety instructions carefully before use, and keep them for future question.

1. Safety indications

- ♥ ALWAYS consult with your physician or supplier to determine proper adjustment and use of the device.
- ♥ Don't store outside and don't let under the sun (Beware : don't let the product in a car in case of high outside temperature !!).
- ♥ Don't use in extreme temperature (above 38°C or 100°F or below 0°C or 32°F).
- ♥ Check carefully all parts of this item.
- ♥ Check regularly if all components are tightened.
- ♥ Use only as a shampoo basin and in proper way.
- ♥ ALWAYS observe the weight limit on the labelling of your product.
- ♥ No modification must be made on this device. Use only original parts to replace the broken or defect one and let this be done by the technical staff.
- ♥ In case you didn't use your product for a long time, let it checked by the technical staff, before using.
- ♥ Do not store/use outdoor, but at home only.
- ♥ Use only on a flat and solid ground.
- ♥ If your skin develops redness, discontinue the use of the cushion immediately and consult your doctor or therapist.
- ♥ Avoid contact with any sharp objects or pointed edges that may

tear or puncture the pad.

2. Installation

Only specialists can adjust the product and teach you how to use it. Ask the technical staff how to use this product and make sure that you understood clearly their explanation.
In case of doubt, don't hesitate to ask again.

Description

The shampoo kit provides an easy, comfortable way to wash the hair of bedridden patients without getting water on the patient or on the linens. The shampoo kit is delivered with an inflatable basin [A] with a drainage hole [E], a drainage hose [E] and a bag [B] to store water to be hang over the basin. Water can be dispensed through a portable showerhead [C] with an on/off valve [D] to control the flow.

Beware : Place the shampoo basin on a flat and semi-rigid surface for better use.

Using

1. Inflate the shampoo basin either by mouth or by using an air pump to get a basin that fits around the client's neck and shoulders comfortably. The air valve has a one-way valve inside to help to prevent air from leaking out. To inflate, squeeze air valve by its neck to help to open the one-way valve inside.
2. Place the drainage hose to a water bucket which must be placed lower than the shampoo basin. Do not punch the hose.
3. Unplug, if necessary, the drainage plug for water drainage while in use.
4. If necessary, hang the water bag on a bedpost or IV pole and use the portable shower to wash hairs. The on/of valve helps controlling the flow.
5. After use, deflate for easy storage.

3. Cleaning

The basin has to be cleaned regularly and after each use with a mild cleaning agent and a dry cloth.

4. Materials

The product is made of sturdy vinyl.

5. Disposal

Dispose of this product in accordance with your regional waste disposal regulations.

6. Technical specifications

This product has a 2 years warranty, from the date of purchase, against any manufacturing defect. This product is biocompatible.

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