GR EN

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INTENDED USE

These wheelchairs have been designed for physically or mentally handicapped people with limited mobility.

It makes possible to transport people with mobility problems or limited mobility and can be used for movement by the patient himself (if it has 20-24" rear wheels) or with the help of an attendant.

When the wheelchair has a rear wheel with a diameter of 50 to 60 cm and hand rim on the wheel then the user can move independently

The Wheelchair with the Commode can also be used as a toilet seat.

SPECIFICATIONS

FRAME	Steel – foldable
FRAME	Single Crossbar
SEAT	Fabric - Nylon
ARM REST	Flip back
FOOTREST	Detachable, Rotating, Swing away with adjustable height
HANDLES	With Adjustable Height
REAR WHEELS	22" MAG wheel with solid tire & Quick release
FRONT WHEELS	6" x 1.25" Solid (150x32mm)
SAFETY BELT	Yes

DIMENSIONS

											
0815009	36cm	35cm	36cm	34 cm	57cm	94cm	87-96* cm	44cm	22"6"	68kg	16**kg

*With the handles in the lowest position, the total height is 87cm. With the handles in the highest position, the total height is 96cm.

**With footrests & rear wheels. The wheelchair without footrests and rear wheels weighs 12.5kg.

ASSEMBLY

Please refer on **ANNEX A**

SYMBOLS

Please refer on **ANNEX B**

GENERAL WARNINGS:

- If you don't read this user manual, it is preferable not to use this product or another available part. If you don't understand the usage or the precautions, please contact with the dealer or the proper technical person before you use the product because may be caused damage.
- Keep away the product from heat sources.
- Be careful when children are nearby and do not allow the children to play with the product.
- Do not exceed the maximum payload.
- Special care is always required where there are moving parts that can cause limb entrapment and injury.
- Do not attempt to lift the wheelchair by any removable parts.
- Getting in and out of the wheelchair may be performed independently, provided you have an adequate mobility and upper body strength.

IMPORTANT!

The construction of this wheelchair is designed to provide maximum safety. However, users may be at risk by using their wheelchairs unsuitable. For your own safety, the following rules must be observed. Nonspecific or incorrect adjustment work can increase the risk of accidents.

- Always wear light-colored clothes when moving in dark places so that you can be reached more easily from others. Be careful during your first course in this wheelchair.
- Meet your wheelchair.
- Make use of your wheelchair on level ground as well as uphill and downhill.
- Do not attempt to reach objects if you need to lift them from the floor reaching between your knees.
- Do not attempt to reach objects if you need to move on to the seat.
- Do not use on roads or motorways.
- Do not attempt to climb or descend with water, ice, or oil membrane.
- Do not use unauthorized parts, accessories, or adapters other than those approved by MOBIAK.
- Do not stand in the wheelchair.

BEFORE USE:

- Check if the frame is damaged in order to ensure the safe use of the product. There are no cracks or breaks in the frame.
- Check if the product is assembled correctly.
- Check that the seat is fully open. If it is not, push the seat down from the sides of the frame.
- Check that the brakes are working properly.
- Always check the state of wear of the mechanical parts in order to ensure the completely safe use of the product for people and objects.
- If there are quick-release wheels, check if they are correctly engaged.

WARNING!

The following methods must be used only from experienced patients. If you are not experienced, let you help from an attendant!

GETTING IN THE WHEELCHAIR:

1. Position the wheelchair as close as possible to you.
2. Engage wheel locks.
3. Lean on the frame with arms and slowly go down towards the seat.

GETTING OUT OF THE WHEELCHAIR:

1. Position the wheelchair as close as possible alongside the seat to which you are transferring.
2. Engage wheel locks.
3. Must rely on both handles equally. Supporting only one hand (and unattended) can cause injury.
4. Shift body weight into seat with transfer.

OTHER INSTRUCTIONS:

1. Always use lifts and ramps. If these are not available, you can reach your destination with the help of two escorts. They must lift the wheelchair by holding stable parts of the frame. A wheelchair should never be lifted while a passenger is seated. It should be pushed only.
2. Avoid stairs or obstacles, otherwise ask help from your attendant:
3. This wheelchair must be used solely to carry a person at a time on the seat.
4. Users should not sit in wheelchair while are using any means of transport.
5. The wheelchairs have a brake on the back wheels, which should be secured when the wheelchair is not moving.
6. DO NOT push on the footrests to sit or stand up, DANGER OF OVERTHROW. The footrests are used by the user when He or She is seated.
7. DO NOT touch moving parts while moving DANGER OF INJURY

MAINTENANCE:

Periodic review of the procedures below will help to maintain the life and effectiveness of your wheelchair.

WEEKLY

1. Check the wheel pressure by following the instructions on the wheels.
2. Check that the central bolt of the transverse arm is tight and in good condition.
3. Clean the fabric with soapy water or sufficient products.
4. Clean the metal parts of the chair with a dry cloth or with special products for chrome surfaces or painted surfaces.
5. Check the wheel bolts and adjust if necessary.
6. Check for damage to front and rear wheels.

QUARTERLY

1. Check and adjust if necessary: wheel forks, pulling wheels, propellers.
2. Grease or oil all components that are subject to friction or impending wear.
3. Tighten the screws and nuts.
4. General revision of the chair's status and repair to whatever is necessary.

ANNUALLY

1. We recommend a full review once a year to a qualified service provider.
2. For the safety of the patient, we recommend that the wheelchair will be checked by the manufacturer every 2 years. In case of repair, use only authentic spare parts.

CLEANING:

In addition to maintenance and inspection, both the frame and the upholstery can be thoroughly cleaned.

- To maintain its appearance, the fabric must be cleaned regularly.
- Dirt can be removed with a sponge and a mild detergent with warm water, sweeping in the direction of texture.
- Rinse with clean water before letting it dry.
- **DO NOT leave the wheelchair wet.**
- Do not place the fabric in a tumble dryer. If the dirt is deep, the surface should be cleaned with a soft brush.
- The frame should be cleaned regularly with a mild detergent.
- The wheelchairs with leather seats must be cleaned with dry cloth.
- Clean the wheels with warm water and mild soap. Pay attention in drying the wheels. it is dangerous use the wheelchair with wet wheels.
- If you need to disinfect the product, use a common, mild disinfectant.
- Protect the item from scratches, cuts and punctures.

CAUTION (!):

- In case of damage or malfunction of your product DOT NOT use the product and please contact the authorized and trained representative.
- Do not use the product for a purpose not indicated in this manual
- MOBIAS S.A declines all responsibilities for any consequences resulting from an incorrect use of this product and from unauthorized alteration to the frame of the product.

- MOBIAC S.A reserves the right to change the information contained in this document without previous notice.
- There is a $\pm 3\%$. Deviation for the Dimensions of the Products.

ENVIRONMENTAL PROTECTION

If one day you find that your product needs to be replaced or it is no longer working for you, consider protecting the environment:

- 1) Do not dispose your product along with the rest of the public waste (this is also the meaning of the shown recycling sign).
- 2) Contact your Public Authorities and they will instruct you of the Recycling centers to which your product must be disposed.
- 3) Correct disposal of your product helps the protection of the environment as well as the recycling of the product's components.

DECLARATION OF CONFORMITY:

We are solely responsible for declaring that the Medical Devices mentioned in this statement are of Low-Risk Class (Class I) and comply with the requirements of the European Regulation 2017/745 and where appropriate, the standards and legislation referred to.

INCIDENT REPORTING:

Any serious incident that occurs in relation to use of the device must be reported to the manufacturer and a competent authority of the Member State in which the user and / or patient resides.

WARRANTY:

The product is guaranteed for a period of two (2) years from the purchasing date.

The warranty covers factory defects & DOES NOT cover damages resulting from misuse, poor maintenance, modification, overuse or non-compliance with the product's operating instructions. Also DOES NOT cover the parts that can wear out during use or over time, like the upholstery, wheels or brakes.

Damages or defects caused by natural disasters, unauthorized maintenance or repair, supply problems (where applicable) & carriers are not covered by the warranty.

No service or spare part is covered unless it has been inspected by an authorized service department or the manufacturer of the product. Shipping costs for goods of warranty are to be paid by the user- customer or by dealer.

WARRANTY CLAIM FORM

BUYERS DETAILS:			
FULL NAME:			
PURCHASING DATE:		PHONE:	
SERIAL NUMBER:			
DEALER DETAILS			
FULL NAME:			
ADDRESS:			
PHONE:			
		STAMP AND SIGNATURE	

MOBIAC S.A
KAGIANA AKPOTHPIOY, XANIA - KPTH
KATHIANA AKROTIRIOU CHANIA - CRETE
T: +30 28210 63222 | WEB: www.mobiackcare.com

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The wheelchair that you have in your hands may differ from the one in the following photos. The assembly process is similar.

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First, we check the contents of the package.

Place the wheelchair as shown in the picture.

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By pushing the bars on the side of the seat with one movement, the wheelchair opens and comes into wide open state.

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
Hold the footrest with your hand and place it in the way shown. The pin must pass through its hole

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Use the included belt to make sure that the patient's feet are properly positioned.

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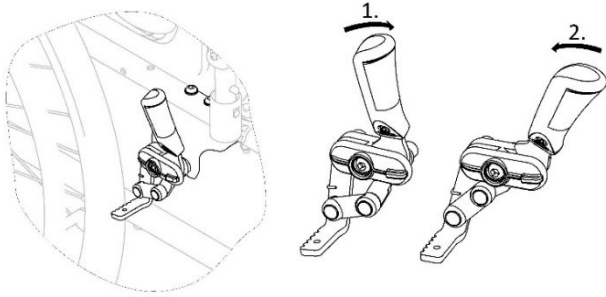The brakes are used to immobilize the wheelchair when it is stationary to prevent it from rolling away.

To engage the brake, push the brake lever forwards with the palm. (1)

To disengage the brake, pull the brake lever backwards. (2)

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Never engage the parking brakes while you are moving. Do not engage the parking brakes to slow down the wheelchair.

Do not lean on the parking brakes for support or transfer. Keep your fingers away from movable parts when using the brake, always keep your hand on the brake lever.

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How to remove the Quick Release rear wheels:

The Wheelchair has a quick-release system on the rear wheels so they can be removed easy and can reduce also the storage volume in the car and home.

To disassemble the wheel, press the black button in the center of the wheel, to unlock, and at the same time pull the wheel outwards until it comes out completely from its base.

Repeat the same process on the other wheel.

To assemble follow the reverse procedure.

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Description of symbols that have been printed on the label, Master Carton and or User Manual

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